

THERMO SHAKER INCUBATOR



Description :-

TS100/TCS10 Thermo Shaker is a temperature controlled shaker, which is useful for a variety of molecular biology applications. Speed, time and temperature settings are continuously visible on the LCD, simultaneously showing both actual and selected values. Integral over-temperature control ensures long time, safety and sample integrity.

FEATURES :-

- LCD display. It is easy to setup and use.
- Reliable design with inner extra temperature protection.
- Each device is equipped with lifter that is easy to replace different blocks.
- Digital accuracy and convenience.
- Precision wells for uniform thermal transfer.
- Assorted blocks for a wide variety of tube sizes.
- Brushless DC Motor: noiseless long service life maintenance-free.
- Conforms to CE safety standard.
- Accurately control and display time, temperature and speed.
- 12-month warranty.



sales@helixbiosciences.com



www.helixbiosciences.com



+91-11-41613622
1800-890-4669

REGISTERED OFFICE :-

A-7, SOUTH WEST DELHI, NARAINA INDUSTRIAL AREA PHASE 1, NEW DELHI - 110028

MANUFACTURING :-

PLOT NO. 360, IMT INDUSTRIAL ESTATE, SECTOR 68, FARIDABAD, HARYANA - 121004

THERMO SHAKER INCUBATOR



Applications :-

- Maintaining 14°C for ligation reactions.(HBSTS10)
- Maintaining 17°C for storing oocytes.(HBSTS10)
- Storing samples at ice bucket temperatures.
- Storing enzymes or DNA libraries at your work station.
- Transporting samples from the field at exactly the right temperature.
- Heating samples in the field.
- Enzyme reactions and deactivations.
- PCR sample preparation.

Applications :-

MODEL	HBSTS100	HBSTS10
Temp setting range	0°C~100°C	0°C~100°C
Temp control range	RT+5°C~100°C	RT.-20°C~100°C
Time range	1min ~ 99h59min/∞	1min ~ 99h59min/∞
Temp accuracy	≤±0.3°C	≤±0.3°C*



THERMO SHAKER INCUBATOR




Applications :-

MODEL	HBSTS100	HBSTS10
Display accuracy	0.1°C	0.1°C
Temp uniformity	≤±0.3°C*	≤±0.3°C*
Speed range	200~1800rpm	200~1500rpm
Orbit	3mm	3mm
Heating time	≤12min(from 25°C to 100°C)	≤15min(from 25°C to 100°C)*
Cooling time1		≤10min(from 100°C to 25°C)*
Cooling time2		≤20min(from RT. to RT.-20°C)
Multi point running	Yes(Max 5points)	Yes(Max 5points)
Multi cycle	Yes(Max 99times)	Yes(Max 99times)
Auto heating	Yes	Yes
Auto running	Yes	Yes
Auto resume to run	Yes	Yes
Storage at 4°		
Power supply	150W	150W
Voltage	AC220V or AC110V/50-	AC110V~220V/50-60Hz
Fuse	250V 3A Φ5×20	250V 3A Φ5×20
Dimension(mm)	260x195x150	260x195x150
Weight	7.5kg	7.0kg



THERMO SHAKER INCUBATOR

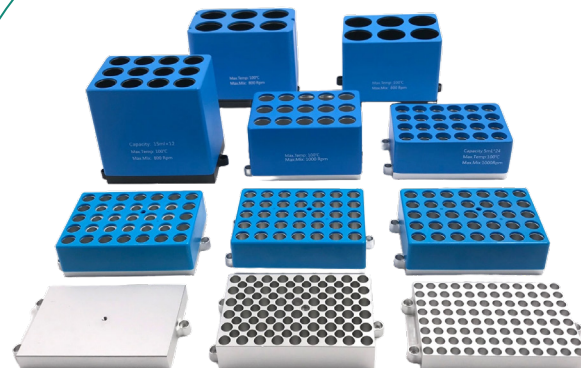

Optional Blocks :-

Type	Capacity	Temp range	Max Speed (HBSTS100)	Max Speed (HBSTS10)
MD30	96x0.2ml	0°~100°C	1800rpm	1500rpm
MD31	54x0.5ml	0°C~100°C	1800rpm	1500rpm
MD32	35x1.5ml	0°C~100°C	1800rpm	1500rpm
MD33	35x2.0ml	0°C~100°C	1800rpm	1500rpm
MD34	20x1.5ml+15x0.5ml	0°C~100°C	1800rpm	1500rpm
MD35	96x0.2ml	0°C~100°C	1800rpm	1500rpm
MD36	24x5ml	0°C~100°C	1000rpm	1000rpm
MD37	15x10ml	0°C~100°C	1000rpm	1000rpm
MD38	12x15ml	0°C~100°C	1000rpm	800rpm
MD39	6x50ml	0°C~100°C	800rpm	800rpm

